


CART FOR DEVICES

BTL-5000 AND BTL-08 ECG

ASSEMBLY INSTRUCTIONS


1 ASSEMBLY INSTRUCTIONS

1.1 LIST OF COMPONENTS AND PARTS

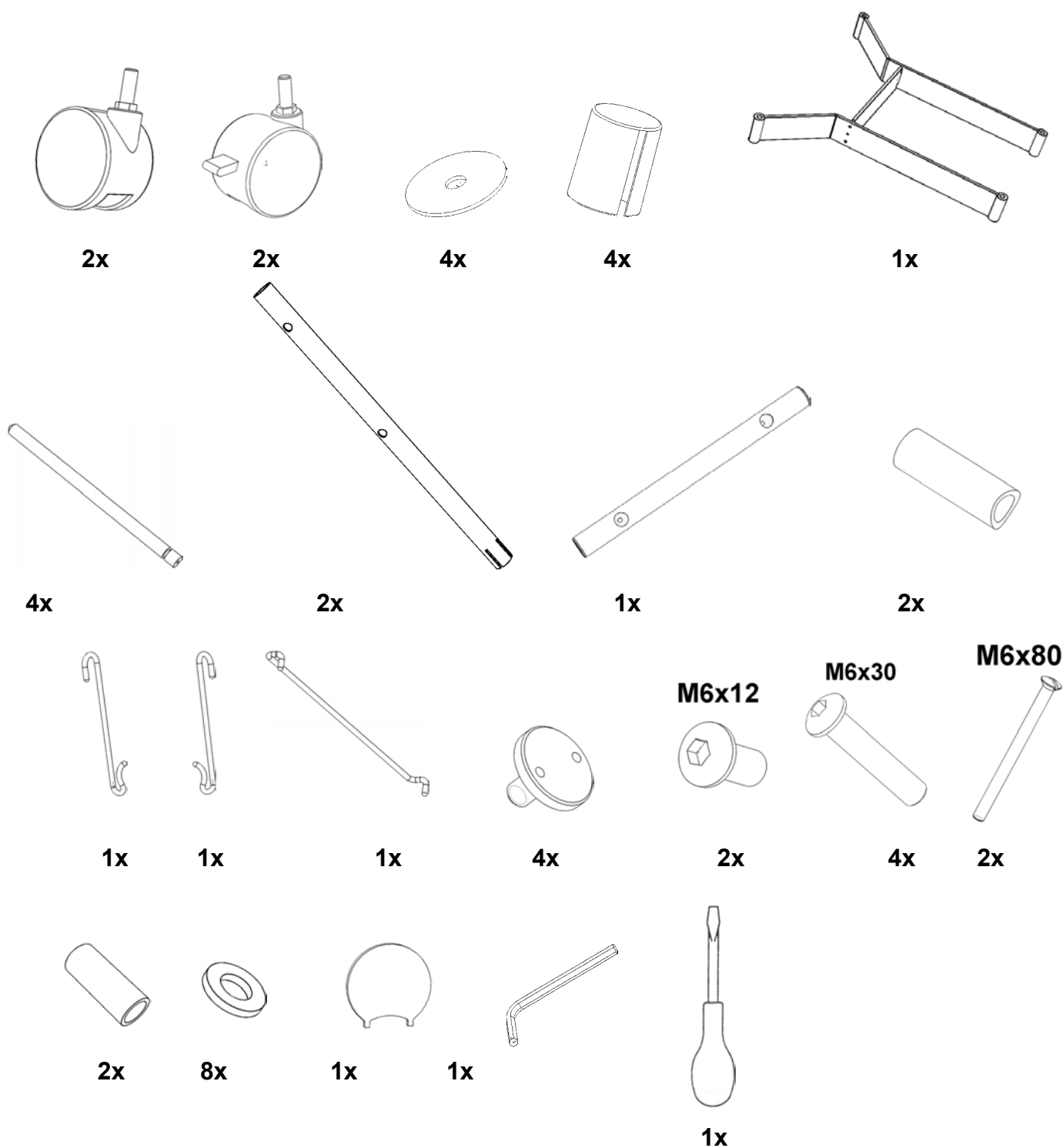

1.2 AVAILABLE MODELS

1.2.1 ECG Cart

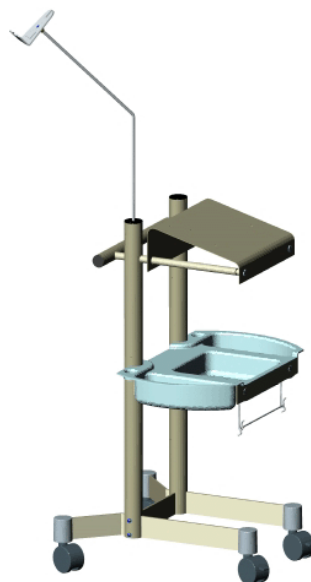


1.2.2 Physiotherapy Cart




1.3 OPTIONAL ACCESSORIES

1.3.1 Attachment for BTL-08 ECG

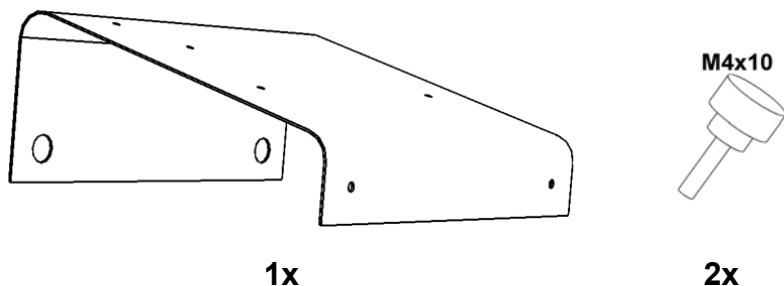


Two hand bolts for fixing the BTL-08 ECG device are delivered together with the ECG attachment.

1.3.2 Attachment for BTL-5000

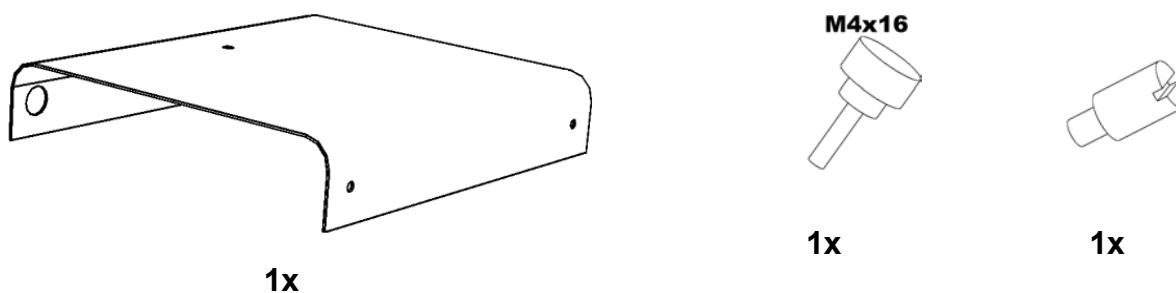


Components for fixing the BTL-5000 unit (1x hand bolt M4x16 and 1x peg) are delivered together with the physio attachment.

1.3.3 Shelf

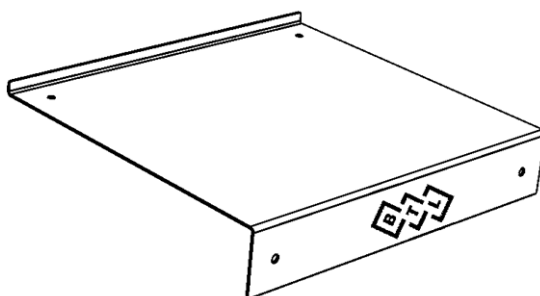


The shelf can be used as storage space or for placing of BTL VAC.

1.3.4 Dish

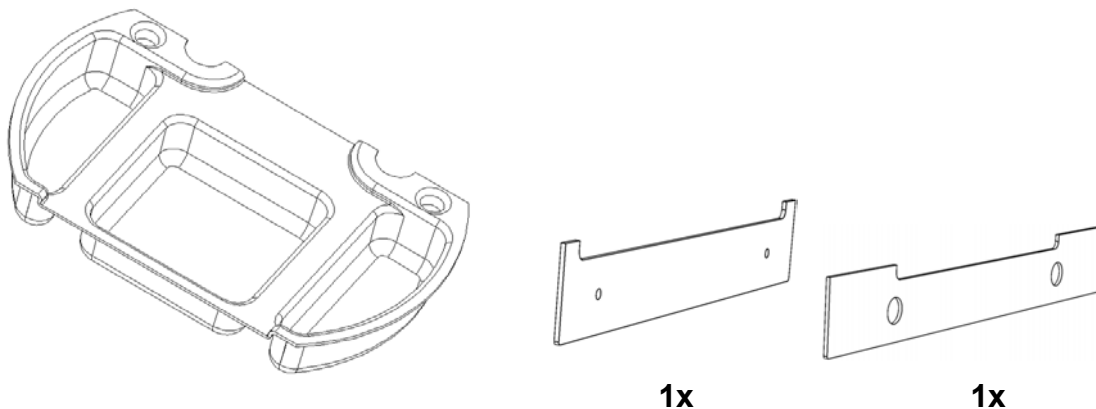



The dish can be used as three-part storage space or as space for placing of BTL VAC with two storage areas at the sides of the BTL Vac device.

1.3.5 Holder for ECG Cable

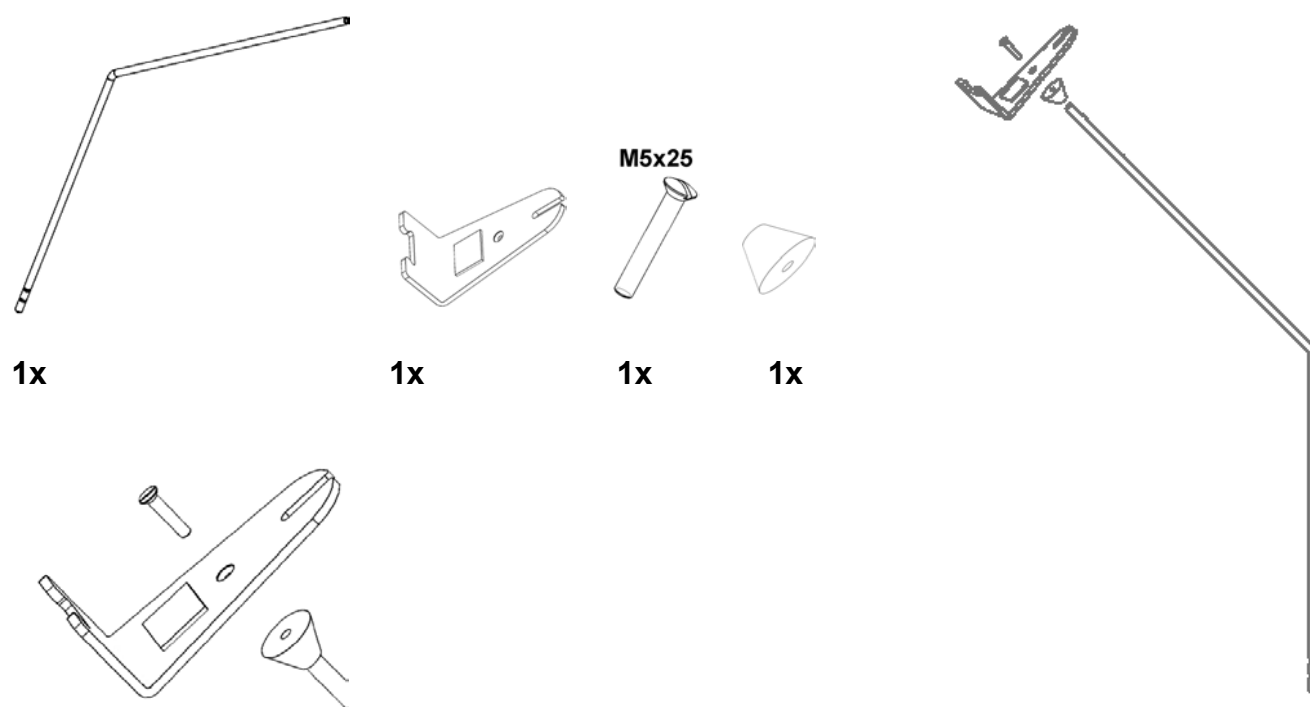


The holder for ECG cable is standardly delivered with the ECG attachment.

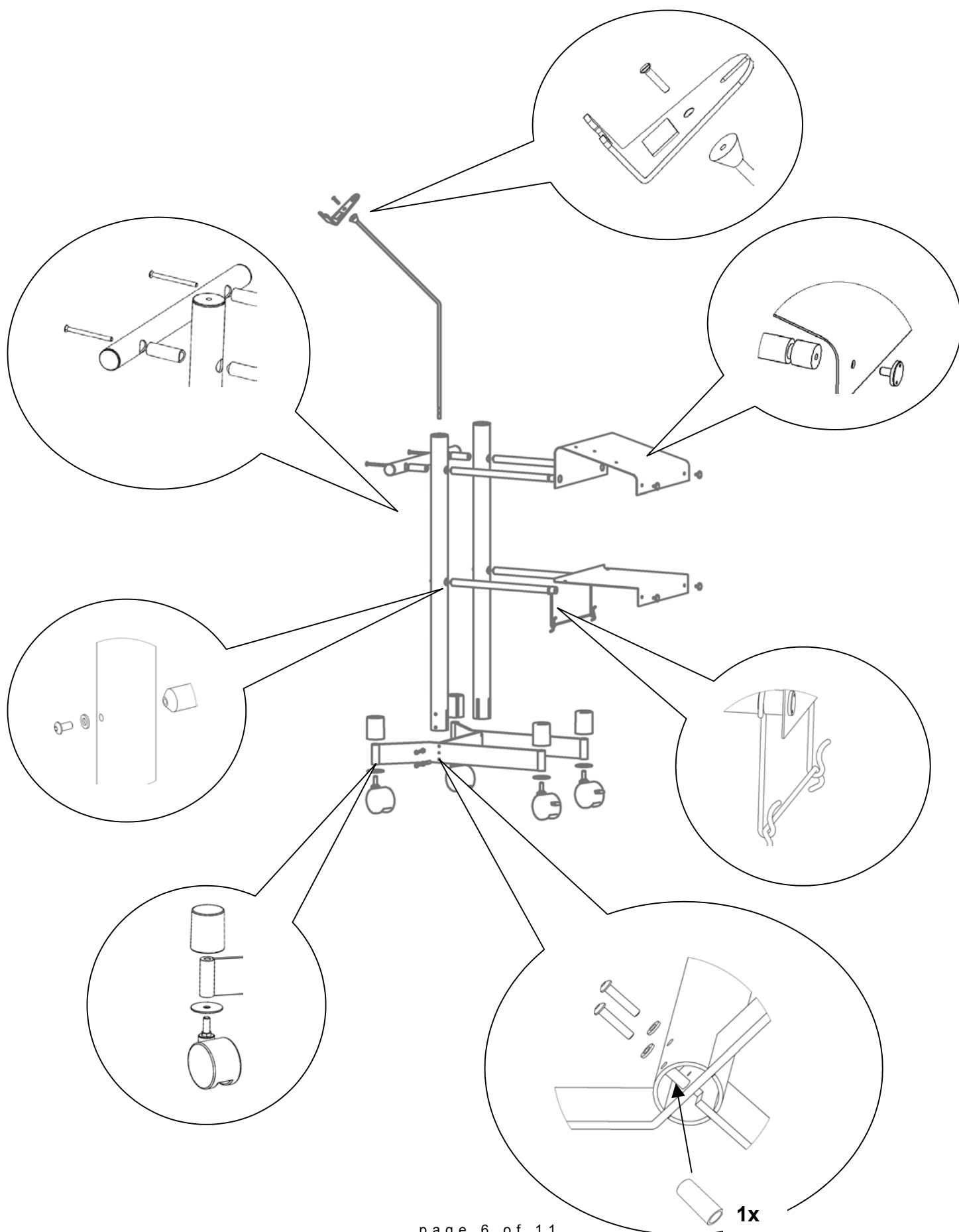

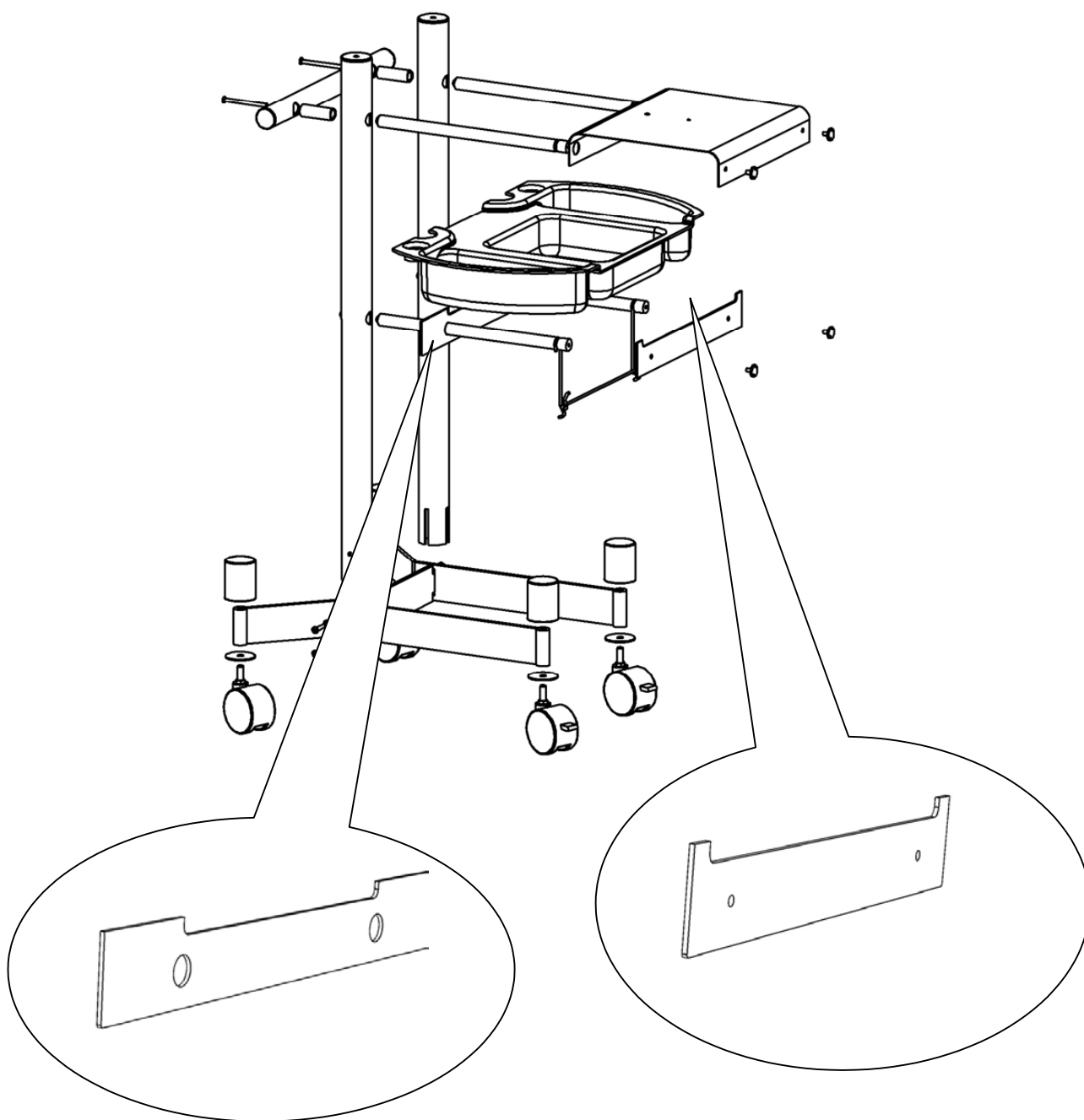


1.4 MOUNTING OF THE BTL-5000 DEVICE TO THE PHYSIO ATTACHMENT

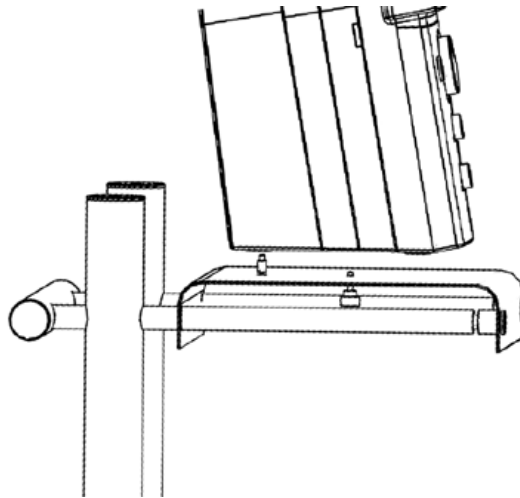

The BTL-5000 devices shall be fixed to the attachment by the M4x16 hand bolt and the peg.

- By the screwdriver screw the peg to the bottom side of the BTL-5000 device.
- Place the BTL-5000 device on the attachment so that the peg fits in the hole in the attachment.
- Straighten the BTL-5000 device so that its front side is parallel with the front edge of the attachment.
- From the bottom of the attachment push the M4x16 hand bolt through the hole and screw it to the device.
- Connect the BTL-5000 device according to the User's Manual.



1.5 MOUNTING OF THE BTL-08 ECG DEVICES TO THE ATTACHMENT

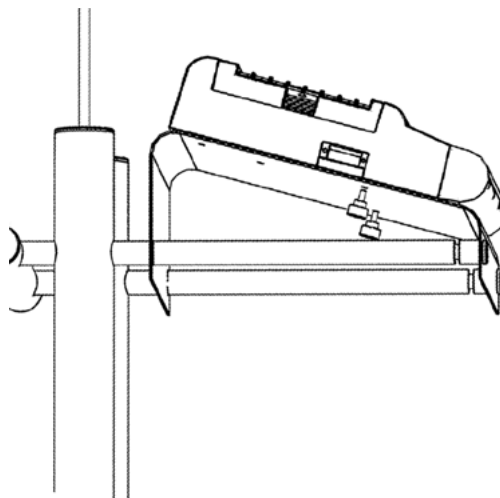

The BTL-08 ECG devices shall be fixed to the attachment by two hand bolts M4x10.

- Place the BTL -08 ECG device on the attachment so that the holes in the attachment match with the thread inserts on the bottom of the device.
- From the bottom of the attachment push the M4x10 hand bolts through the holes and screw them to the device.
- Connect the BTL-08 ECG device according to the User's Manual.

1.6 PLACING OF BTL -VAC

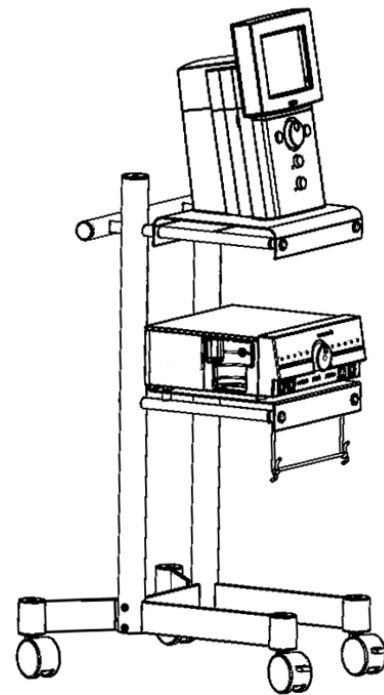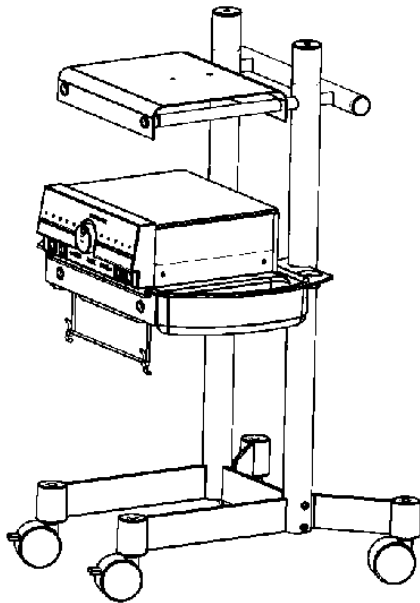

The BTL Vac device shall be placed on the shelf or dish in its clearly outlined place.

- Unpack the BTL-Vac device.
- Place it on the shelf or dish.
- The front bottom edge of BTL-Vac shall be aligned with the front edge of the shelf or plastic dish.
- Connect BTL-Vac according to the User's Manual.

1.7 MAINTENANCE

All movable parts of the device are manufactured and finished in the way that in case of standard use it is not necessary to maintain (grease) them during the whole lifetime of the cart.

1.7.1 General Principles of Cleaning

Foreign substances (paint, etc.) on the surface shall be removed by a cloth and tepid water. Never use **steel wool, metallic sponges, knives or scouring washcloth**. Stains of oil or grease can be removed by **denatured alcohol**.

Never use strong solvents such as **acetone, paint thinners, benzine, ammonia or chloride agents**, or **abrasive polish** – otherwise the finishing or plexiglass could be damaged.

1.8 PRINCIPLES OF SAFE HANDLING

- Do not place other devices on the cart! The cart is designed exclusively for BTL devices. Placing of other devices or equipment can cause their damage or endanger the personnel and/or the patient!
- The devices to be used together with the cart (BTL-08 ECG and BTL-5000) must be always securely fixed to the attachment, so that they cannot get loose and fall.
- When moving the cart make sure that all conductors of the devices are securely attached and will not be dragged on the ground. The cables of electrodes of BTL-Vac shall be looped round the handle of the cart.
- When moving the cart the castors must be unlocked.
- After placing the cart on the final place lock the castors so as to prevent spontaneous or undesired moving of the cart.
- Do not tilt the cart, otherwise the components on the cart or in storage spaces can get loose, fall and suffer damage!
- For instructions for use, secure handling, maintenance and cleaning of BTL devices see the User's Manuals which are supplied together with the devices.

Manufacturer

Medical Technologies, s.r.o.
Heleny Malirove 11
169 00 Praha 6

Sales Department

BTL, s.r.o.
Heleny Malirove 11
169 00 Praha 6
Czech Republic
Tel.: (+420) 220 517 012, 220 514 753
Fax: (+420) 220 515 520
E-mail: sales@btlnet.com
<http://www.btlnet.com>

© All rights reserved. No part of this manual may be reproduced, saved in a research centre or transferred by any means incl. electronic, mechanic, photographic or other records without previous approval from MEDICAL TECHNOLOGIES, s.r.o.

BTL is a registered trademark of MEDICAL TECHNOLOGIES, s.r.o.

MEDICAL TECHNOLOGIES, s.r.o operates a policy of continuous development. Therefore, it reserves the right to make changes and improvements to the Product described in this manual without prior notice.

The contents of this document are provided "as is". Except as required by applicable law, no warranties of any kind, either expressed or implied, are made in relation to the accuracy, reliability or contents of this document. MEDICAL TECHNOLOGIES, s.r.o reserves the right to revise this document or withdraw it at any time without prior notice.